

TransIT[®]-AAViator Transfection System



Protocol for MIR 73750 and MIR 73745

Quick Reference Protocol, SDS and Certificate of Analysis available at mirusbio.com/literature

INTRODUCTION

Adeno-associated virus (AAV) is a nonenveloped, single-stranded DNA virus from the *Parvoviridae* family notable for its lack of pathogenicity, low immunogenicity and ability to infect both dividing and quiescent cells. Because AAV is replication-defective in the absence of adenovirus or helper proteins and is not implicated in any known human diseases, it is widely considered a safe gene delivery vehicle for *in vivo* and *in vitro* applications. Accordingly, recombinant AAV has become an invaluable tool for gene therapy and the creation of isogenic human disease models.

The TransIT[®]-AAViator Transfection Reagent enables the generation of high titer AAV in HEK 293 cell types while improving the full/empty capsid ratio and reducing the amount of plasmid DNA required per transfection. The TransIT[®]-AAViator Transfection System includes RevIT[™] AAV Enhancer to further enhance the performance of TransIT[®]-AAViator Transfection Reagent, making this kit ideal for generating high-titer AAV preparations to accelerate research and development.



For best results, TransIT[®]-AAViator Transfection Reagent is intended to be used with RevIT[™] AAV Enhancer.

SPECIFICATIONS

Storage	Store TransIT [®] -AAViator Transfection Reagent and RevIT [™] AAV Enhancer at -10 to -30°C, tightly capped. <i>Before each use</i> , warm to room temperature and vortex gently. RevIT [™] AAV Enhancer is known to maintain function through at least five freeze-thaw cycles (thawed in a 37°C incubator). Return to proper storage conditions after each use.
Stability / Guarantee	When properly stored and handled, TransIT [®] -AAViator Transfection Reagent and RevIT [™] AAV Enhancer are guaranteed for 6 months from date of purchase.



Warm all reagents to room temperature and mix gently before each use.

NOTE: RevIT[™] AAV Enhancer remains frozen at temperatures < 19°C.

MATERIALS

Materials Supplied

The TransIT[®]-AAViator Transfection System is supplied in the following formats:

Product No.	Volume of TransIT [®] -AAViator Transfection Reagent	Volume of RevIT [™] AAV Enhancer
MIR 73750	2 × 1.5 ml	1 × 1.5 ml
MIR 73745	1 × 30 ml	10 × 1.5 ml

For Materials Required, but Not Supplied, See Page 3.

For Research Use Only

BEFORE YOU START:

Important Tips for Optimal AAV Production

The suggestions below yield high-efficiency plasmid DNA transfection using the *TransIT®-AAViator* Transfection System.

- **Cell culture conditions.** Culture suspension HEK 293 cells in appropriate complete growth medium. Before transfection, ensure cells are $\geq 95\%$ viable by trypan blue exclusion (or similar method) and doubling every 24 hours. After transfection, there is no need to perform a medium change to remove the transfection complexes.
- **Cell density at transfection.** The recommended cell density for suspension cells is 3×10^6 cells/ml. Passage cells 18-24 hours before transfection to ensure that cells are actively dividing and reach the appropriate density at time of transfection.
- **AAV packaging and transfer plasmids.** The optimal ratio between plasmids will depend on the vector backbone and gene-of-interest. For each unique construct, empirically determine and use the optimal ratio for best results. Use plasmid manufacturer recommendations or previously established ratios as a starting point.
- **Ratio of *TransIT®-AAViator* to DNA.** The optimal ratio between *TransIT®-AAViator* Reagent and DNA is typically 1:1 to 1.5:1, therefore, the recommendation is to start at 1.25:1. Determine the optimal *TransIT®-AAViator* Reagent:DNA ratio for each cell type by varying the amount of reagent from 1-2.3 μ l per 1-1.5 μ g of total DNA. Refer to **Table 1** for recommended starting conditions based on culture size.
- **RevIT™ AAV Enhancer.** Titrate RevIT™ AAV Enhancer from 0.5 to 1.5 μ l per 1 ml of culture to determine the optimal amount for production of your specific viral vector.
- **Complex formation conditions.** Prepare *TransIT®-AAViator*:RevIT™ AAV Enhancer:DNA complexes in compatible basal cell culture media in a volume that is 5% of the total culture volume. We recommend a complex formation time of 15-45 minutes. If forming complexes in a volume that is greater or less than 5% of the total culture volume, complex formation time may need adjustments.




Premix packaging and transfer plasmids together prior to adding to the complex formation medium.



Do not use serum or antibiotics in the media during transfection complex formation.

Transfection complexes can be added directly to cells cultured in growth media +/- serum and up to 0.1-1X antibiotics.

TransIT®-AAViator Transfection System Workflow:

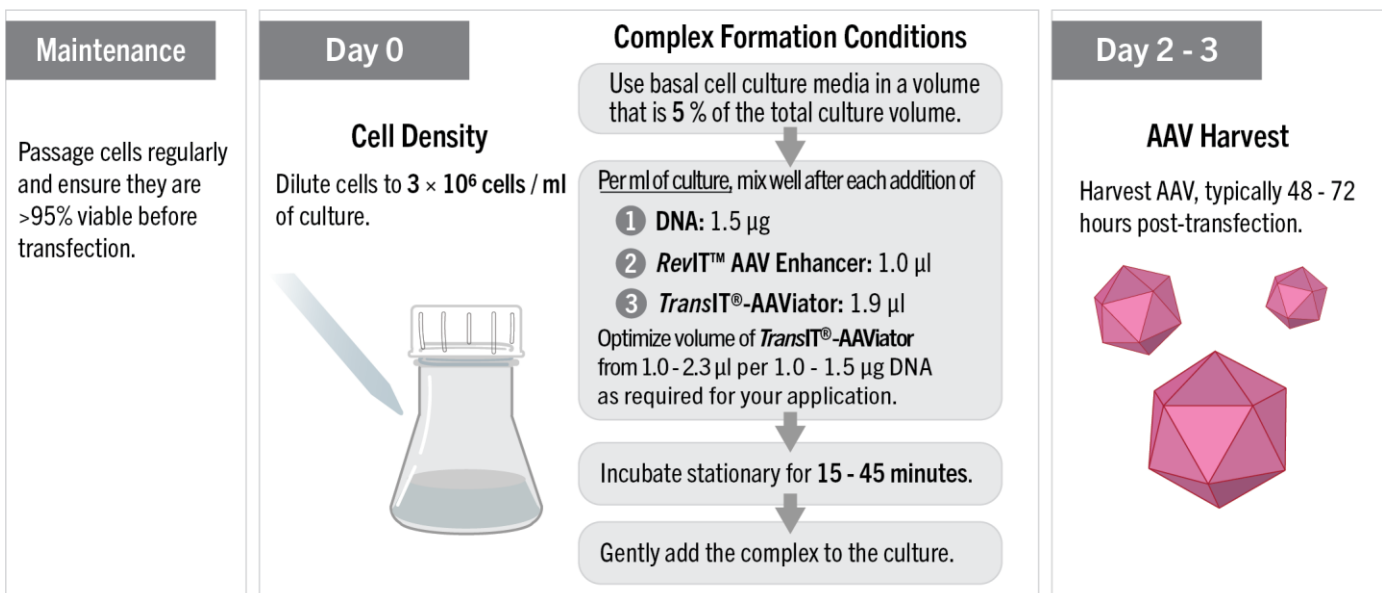

AAV GENERATION IN SUSPENSION HEK 293 CELL CULTURES

NOTE: Use of the *TransIT®-AAViator* Transfection System is only recommended for AAV production in suspension HEK 293 cell lines. Contact Mirus Bio Technical Support for optimization in adherent cell culture platforms.

The following procedure describes plasmid DNA transfections for AAV generation in 125 ml Erlenmeyer shake flasks using 30 ml of complete growth medium. If using an alternate cell culture vessel, increase or decrease the amounts of *TransIT®-AAViator*, *RevIT™* AAV Enhancer and total DNA based on the **volume of complete growth medium** to be used. To calculate the required reagent quantities based on the recommended starting conditions and total culture volume, refer to the calculation worksheet in **Table 1** (below).

Table 1. Scaling worksheet for *TransIT®-AAViator* with *RevIT™* AAV Enhancer

Starting conditions per milliliter of complete growth medium				
	Per 1 ml		Total culture volume	Reagent quantities
Basal Culture Medium	0.05 ml	×	_____ ml	= _____ ml
Total Plasmid DNA (1 µg/µl stock)	1.5 µl	×	_____ ml	= _____ µl
<i>RevIT™</i> AAV Enhancer	1.0 µl	×	_____ ml	= _____ µl
<i>TransIT®-AAViator</i> Reagent	1.9 µl	×	_____ ml	= _____ µl

NOTE: Total Plasmid DNA refers to the combined weight of AAV plasmids (in µg) per transfection. For example, pre-mix AAV packaging and transfer plasmids to a stock solution of 1 µg/µl and add 1.5 µl of this mixture to deliver 1.5 µg of Total Plasmid DNA.

Materials Required, But Not Supplied

- Suspension HEK 293 Cells (e.g. Viral Production Cells 2.0, Gibco Cat. No. A49784)
- Complete Culture Medium (e.g. Viral Production Medium, Gibco Cat. No. A4817901)
- Plasmid DNA (e.g. pALD-ITR-WPRE-GFP (Aldevron Cat. No. 5069-10), pALD-HELP (Aldevron Cat. No. 5082-10), AAV8 Rep-Cap Plasmid (GeneMedi Cat. No. P-RC09))
- Erlenmeyer shake flasks (e.g. Corning® Cat. No. 431143 or Thomson Cat. No. 931110)
- 10X Cell Lysis Buffer (500 mM Tris pH 8, 10% Tween® 20, 20 mM MgCl₂)
- 5 M Sodium Chloride (5 M NaCl)
- Benzonase® or equivalent (e.g. Sigma Cat. No. E1014 or Syd Labs Cat. No. BP4200)

Transient Transfection Protocol per 30 ml HEK 293 Culture

A. Maintenance of cells

1. Passage suspension HEK 293 cells 18-24 hours prior to transfection to ensure that cells are actively dividing at the time of transfection and to obtain a density of $3 - 4 \times 10^6$ cells/ml the next day.
NOTE: Perform cell counts and evaluate viability daily to ensure that cells are doubling every 24 hours and $\geq 95\%$ viable by trypan blue exclusion. Do NOT proceed with transfection if cells are not doubling normally or are $< 95\%$ viable.
2. Incubate cells overnight under appropriate conditions (e.g. 37°C , 5-8% CO_2 , shaking).

B. Prepare TransIT®-AAViator:RevIT™ AAV Enhancer:DNA complexes (immediately before transfection)

1. Warm TransIT®-AAViator Reagent and RevIT™ AAV Enhancer to room temperature and vortex gently before using. RevIT™ AAV Enhancer can be incubated in an incubator that is set to 37°C to accelerate thawing. If thawing at room temperature, allow ~ 4 hours and ensure ambient air temperature is $> 19^\circ\text{C}$.
2. Immediately prior to transfection, seed cells at a density of 3×10^6 cells/ml into a transfection culture vessel (e.g. 30 ml per 125 ml Erlenmeyer shake flask).
3. Place 1.5 ml of basal culture medium in a sterile tube.
4. In a separate sterile tube, combine AAV plasmids per manufacturer recommendations to a final concentration of $1 \mu\text{g}/\mu\text{l}$. Mix thoroughly.
5. Transfer 45 μl of the DNA mixture prepared in Step B.4 to the tube containing basal culture medium. Mix completely.
6. Add 30 μl of RevIT™ AAV Enhancer to the diluted DNA and basal culture medium. Mix completely.
7. Add 57 μl of TransIT®-AAViator Reagent to the diluted DNA:RevIT™ mixture. Mix completely by inversion or vortexing. Do **NOT** agitate TransIT®-AAViator:RevIT™:DNA complexes again after this initial mixing.
8. Incubate at room temperature for 15-45 minutes without additional agitation to allow transfection complexes to form.

C. Distribute the complexes to cells in complete growth medium

1. Add the TransIT®-AAViator:RevIT™:DNA complexes (prepared in Step B) to culture vessel, swirling gently to distribute.
2. Shake flasks on an orbital shaker (125 rpm when using a shaker with a 1.9 cm orbital throw) at appropriate temperature and CO_2 levels (e.g. 37°C , 5-8% CO_2).
3. Incubate cultures for 48-72 hours prior to AAV harvest.

D. Harvest and storage of AAV

1. Following the 48-72 hour incubation, transfer the total volume of cell suspension (i.e. ~ 32 ml) to a sterile conical tube or appropriate vessel.
2. Add 0.1X volume of 10X Cell Lysis Buffer (i.e. 3.2 ml) and 100 U/ml Benzonase® (i.e. 3,200 U). Mix completely and incubate at 37°C for 1.5 hours with shaking.
3. Add 0.1X volume of 5 M NaCl (i.e. 3.2 ml) and mix completely. Incubate at 37°C for 30 minutes with shaking.
4. Centrifuge the mixture at $4,100 \times g$ for 10 minutes to remove cell debris. Carefully transfer the AAV containing supernatant to a new sterile tube.
5. Store AAV stocks at -80°C .



Divide cultured cells 18-24 hours before transfection to ensure that the cells are actively dividing at the time of transfection.



Do NOT allow the TransIT®-AAViator to incubate alone in complex formation solution > 5 minutes, i.e. if the reagent is pre-diluted, add DNA within 5 minutes for optimal complex formation.

Do NOT agitate TransIT®-AAViator:RevIT™:DNA complexes after the initial mixing. This will result in decreased titer.



There is no need to change culture medium after transfection, unless required by your cell type or culture conditions.



Benzonase® is a non-specific endonuclease used to liberate virus particles from residual nucleic acids in the cell lysates and increase AAV titers.

TROUBLESHOOTING GUIDE

POOR DNA TRANSFECTION EFFICIENCY	
Problem	Solution
Incorrect vector sequence	If you do not observe expression of your target insert, verify the sequence of your plasmid DNA.
Suboptimal <i>TransIT</i> ® Reagent: <i>RevIT</i> ™:DNA ratio	Determine the best <i>TransIT</i> ®-AAViator Reagent: <i>RevIT</i> ™ AAV Enhancer:DNA ratio for each cell type. Titrate the <i>TransIT</i> ®-AAViator Reagent volume from 1-1.5 µl per 1 µg of DNA. Titrate the <i>RevIT</i> ™ AAV Enhancer volume from 0.5-1.5 µl per 1 ml of culture. Refer to “Before You Start” on Page 2 for recommended starting conditions.
Suboptimal DNA concentration	Determine the DNA concentration accurately. Use plasmid DNA preps with an A _{260/280} of 1.8-2.0. The optimal DNA concentration generally ranges between 1.0-1.5 µg per 1 ml of culture. Start with 1.5 µg DNA per 1 ml of culture. Consider testing different amounts of DNA while scaling the amount of <i>TransIT</i> ®-AAViator accordingly.
Low-quality plasmid DNA	Use highly purified, sterile, endotoxin- and contaminant-free DNA for transfection.
Cells not actively dividing at the time of transfection	Use cesium chloride gradient or anion exchange purified DNA which contains levels of endotoxin that do not harm most cells.
Time of AAV harvest not optimal	Divide the culture at least 18-24 hours before transfection to ensure that the cells are actively dividing and reach optimal cell density at time of transfection. DO NOT proceed with transfection if cells are not doubling normally or are < 95% viable by trypan blue exclusion.
<i>TransIT</i> -AAViator™ was not mixed properly	Determine the optimal time to harvest AAV post-transfection. Though typically 48-72 hours post-transfection, the best time to harvest will depend on the vector construct and production platform.
Disruption of transfection complex formation	Warm <i>TransIT</i> ®-AAViator Reagent to room temperature and vortex gently before each use.
Precipitate formation during transfection complex formation	If <i>TransIT</i> ®-AAViator Reagent is pre-diluted in complex formation solution, DNA should be added within 5 minutes. Incubating the <i>TransIT</i> ®-AAViator Reagent in complex formation solution alone for an extended time results in reduced production of functional virus.
Proper experimental controls were not included	After initial mixing of DNA, <i>RevIT</i> ™ AAV Enhancer and <i>TransIT</i> ®-AAViator Reagent, do not agitate the Reagent:Enhancer:DNA complexes again, e.g. do not vortex or invert before adding to cultures.
	During complex formation, scale all reagents according to the table in the protocol, including serum-free media, <i>TransIT</i> ®-AAViator Reagent, <i>RevIT</i> ™ AAV Enhancer and plasmid DNA.
	Precipitation may be observed when excess DNA is used during complex formation. This may negatively impact transfection efficiency. To avoid precipitation when using high concentrations of DNA, increase the volume of serum-free medium during complex formation by two-fold.
	Large-volume transfection complexes may appear turbid – typically, this phenomenon does <i>not</i> negatively impact transfection as long as complexes are well mixed.
	To assess delivery efficiency of plasmid DNA, use Mirus <i>Label IT</i> ® Tracker Intracellular Nucleic Acid Localization Kit to label the target plasmid or use Mirus pre-labeled <i>Label IT</i> ® Plasmid Delivery Controls (please refer to “Related Products” on Page 7).
	To verify efficient transfection, use <i>TransIT</i> ®-AAViator Reagent to deliver a positive control such as a luciferase, beta-galactosidase or green fluorescent protein (GFP) encoding plasmid.

TROUBLESHOOTING GUIDE continued

HIGH CELLULAR TOXICITY	
Problem	Solution
Cell density not optimal at time of transfection	High toxicity and cell death may be observed if cells are not dense at the time of transfection. For high virus titers using <i>TransIT®</i> -AAViator Reagent, ensure that cell cultures are approximately 3×10^6 cells/ml (for suspension cell transfections) at the time of transfection.
Cell morphology has changed	When generating AAV with <i>RevIT™</i> AAV Enhancer, cell growth may decrease. This is normal and does not adversely affect virus titers.
	Mycoplasma contamination can alter cell morphology and affect transfection efficiency. Check your cells for mycoplasma contamination. Use a fresh frozen stock of cells or use appropriate antibiotics to eliminate mycoplasma.
	A high or low cell passage number can make cells more sensitive and refractory to transfection. Maintain adherent or suspension HEK 293 cells below passage 30 for optimal recombinant virus production.
Transfection complexes not evenly distributed after complex addition to cells	Add transfection complexes while swirling the flask. If this is not possible, gently mix the culture vessel to ensure even distribution of the transfection complexes. However, avoid vigorous agitation that could disturb formed transfection complexes, e.g. vortexing after the initial mixing of the DNA, enhancer and transfection reagent.



Reagent Agent®

Reagent Agent® is an online tool designed to help determine the best solution for nucleic acid delivery based on in-house data, customer feedback and citations.

Learn more at: mirusbio.com/ra

RELATED PRODUCTS

- *TransIT*-VirusGEN® GMP Transfection Reagent
- *TransIT*-VirusGEN® Transfection Reagent
- *RevIT*™ GMP AAV Enhancer
- VirusGEN® GMP AAV Transfection Kit
- VirusGEN® AAV Transfection Kit
- *Label IT*® Plasmid Delivery Controls
- *Label IT*® Tracker™ Intracellular Nucleic Acid Localization Kits
- MiraCLEAN® Endotoxin Removal Kits

For details on the above-mentioned products, visit mirusbio.com.

Contact Mirus Bio for additional information.

Mirus Bio LLC

5602 Research Park Blvd, Ste 210

Madison, WI 53719

Toll-free: 888.530.0801

Direct: 608.441.2852

Fax: 608.441.2849

©1996-2025 All rights reserved. Mirus Bio LLC. All trademarks are the property of their respective owners.

For terms and conditions, visit www.mirusbio.com

Rev1 02Oct2024

Page 7 of 7